

Carodoor

EN

How to operate the Carodoor

When a wireless remote control fails

Operate the Carodoor using the wheelchair lift's original hand control.

In the event of a complete power loss




Disconnect the wheelchair lift from the Carodoor unit.

1. Remove the locking clip attached to the Carodoor link arm by pulling the red loop.
2. Locate the red pry bar (should be attached to the lift pump housing). Place the pry bar between the door rod and motor arm and lever towards the motor housing to release the door rod.
3. Open the doors manually.

If connected to the vehicle's central locking system

1. Unlock the vehicle.
2. You have 4 minutes* to activate the lift before it loses power.
3. If the time has run out and the lift is powerless, unlock the vehicle again or use the lift's wired hand control to activate the system.

*The 4 minute countdown is from version 3.2 of the circuit board. Earlier versions have a 1 minute countdown.

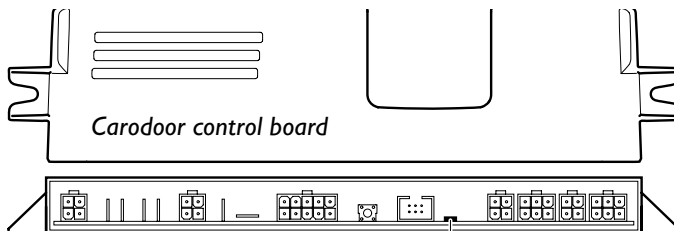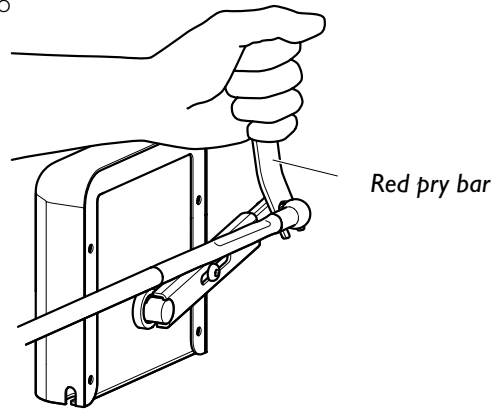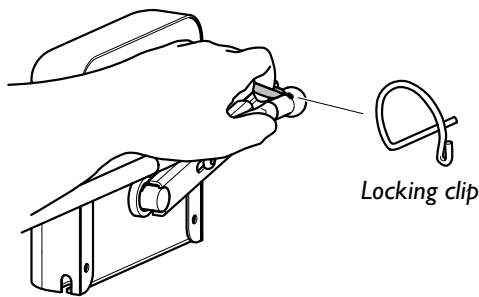

Carodoor status LED-codes.

The Carodoor control board is located inside the lift pump module.

Priority	Color	LED-light	Description
3	Red	0,1 sec (8 pulses)	Lift not secured
2	Red	0,1 sec (1 pulse)	Emergency stop
1	Blue	0,1 sec (1 pulse)	Key pressed
0	Green	0,1 sec (1 pulse)	Lift parked / doors closed
	Green	0,1 sec (8 pulses)	Lift parked / doors in between
	Green	0,1 sec (1 pulse)	Doors open

427484-2 2021/05

Product certificate

Manufacturer /

Braunability Europe AB
Åkerivägen 7, 443 61 Stenkullen
Sweden

Phone: +46 302 254 00
E-mail: info@braunability.eu
www.braunability.eu

Declares that the product /

Carodoor

Conforms to following directives standards and regulations /

2006/42/EEC
ECE R-10

BraunAbility Europe AB, Stenkullen
2019-04-02



Jan M Jensen
Managing Director

EN

Product certificate