

# TILDA™

tilts the seat and makes entering the car much easier



*By keeping the feet on the footrest they won't need lifting.*



*Tilda in combination with the Carony.*

## Tilts the seat in 3 positions

TILDA combines excellently with the seat swivel TURNOUT and makes it even easier to get in and out of the car. By keeping the feet on the footrest they automatically follow over the threshold.

In combination with CARONY, TILDA will provide more head space as the seat slides over from the Carony wheel unit to the swivel base inside the car.

 **auto ADAPT**  
for your independence

# INVENTIONS FOR QUALITY OF LIFE




# TILDA™

TILDA is specially designed to be combined with:  
 TURNOUT swivel base that swivels the seat out over the threshold.  
 and  
 CARONY which enables the seat to be placed on a wheel unit and become a wheelchair.  
 or  
 GLIDE RAILS from Scheelman, providing length adjustability for the seat.

Please note that the footrest is an option.



*Tilda in position 3 and in combination with the TURNOUT swivel base, with the Carony rails mounted on the sides. The footrest is an option.*

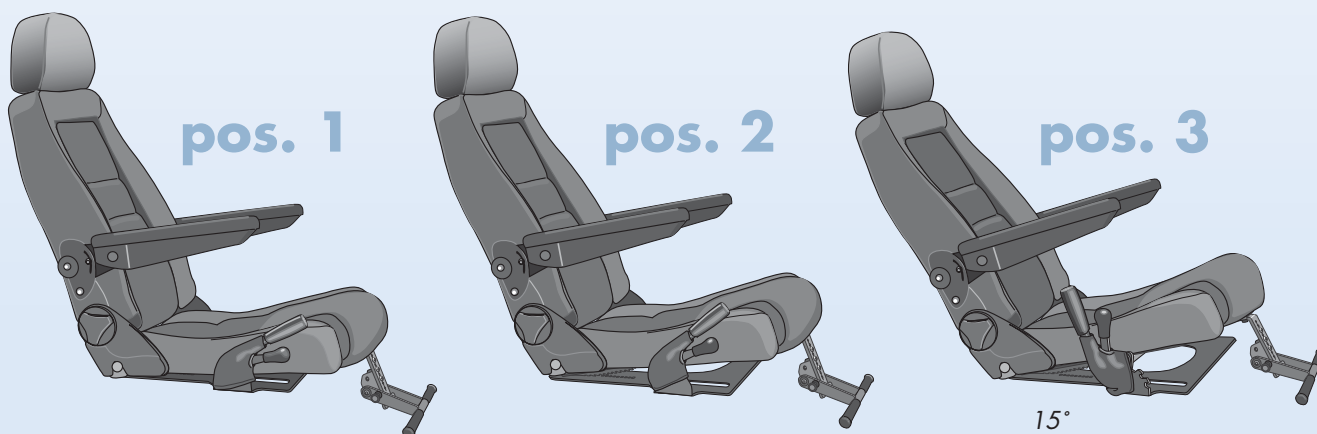

*Grab hold of both levers to release the position of the seat, then raise the seat to an upright position.*

## Crash tested

When tilted in position 1 the ride will become more comfortable and the sitting more relaxed.

TILDA is crash tested in pos. 0 (upright position) and in pos. 1

TILDA is in compliance with following regulations and directives:

74/60/EEC - 2000/4/EG, 74/408/EEC - 2005/39/EEC  
 93/42/EEC, 95/28/EG, ECE R-17, ECE R- 21,  
 SS-EN 12182, ISO 3795  
 FMVSS 201, 207, 208 and 302

Maximum user weight: 150 kg.

Height, under the seat: 24 mm

Width: 520 mm

 **auto  
ADAPT**  
for your independence

**Autoadapt AB**, Åkerivägen 7, S-443 61 Stenkullen, Sweden  
 Phone: +46 302 25400 E-mail: [contact@autoadapt.se](mailto:contact@autoadapt.se)  
[www.autoadapt.com](http://www.autoadapt.com)